



Available on
App store



Available on
Android



iFever thermometer

www.fiismart.com

Product details

Adjustable length	150mm-205mm
Weight	About 20 gms(without battery)
Range	32°C-43°C
Accuracy	$\pm 0.1^{\circ}\text{C}$
Unit	$^{\circ}\text{C}$ or F
Battery	CR2032 3V (210mAh)
Working time	4000hours
Link distance	About 10M(without battery)
Connection mode	BLE4.0
App Software	iFever

Product certification



FCC Certification



CE Certification



RoHS Certification

Safety information

- 1、 Do not drop or twist the thermometer.
- 2、 Do not allow it to be in contact with liquid.
- 3、 The display panel can be disinfected by alcohol, while the part above the display panel can be cleaned by a dry cloth. Do not put it into water or disinfect it in boiling water, it is not waterproof.
- 4、 Do not disassemble the product.
- 5、 The thermometer can be damaged if stored in a high temperature or exposed directly to the sun. Do not keep the product in contact with any chemical solvents in order to avoid dysfunctions caused by chemical reactions.
- 6、 The battery must be recycled or disposed of separately from household waste, to protect human health, water resource and the environment.

Acknowledge

First of all, thank you for choosing iFever smart thermometer. Please read this manual before using the iFever. It's easy for you to operate. The main function is 24h real-time monitoring temperature and high temperature alarm.


APP Introduction

- Thermometric
- High temperature alarm
- Medical records
- Temperature curve
- Cloud backup





- Passed fabric safety inspection

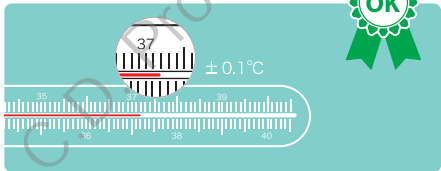

- Passed measurement accuracy test

Warm notification

1. If the received device has any manufacturing defects it must be sent freight prepaid and when the package arrives back to us, we will give a refund which also includes the return shipping costs.
2. The returned package should come along with a copy of the original invoice and the return card with return reasons and the number of returned goods.
3. iFever's warranty does not cover product malfunction arising from:
 - (1) The warranty expiry;
 - (2) Damages resulting from operation & maintenance which are not according to the user guide;
 - (3) Damages caused by misuse, fire damage or immersion in water.
 - (4) Abrasion, corrosion and discoloration after regular use.

Note: This product uses wireless Bluetooth to communicate with cell phones, which is not always stable and reliable under all circumstances. Also, the app data from the related iFever App cannot replace clinical diagnostics.

Operation manual

- 1 Download app 2 Switch On device 3 Wear



iOS



Android




- 4 Open Bluetooth 5 Connect in APP



- 6 Start



Function keys

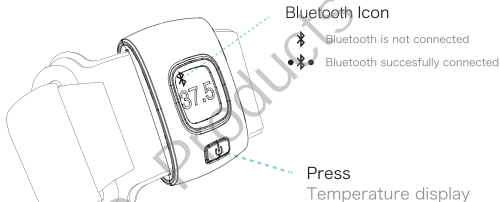

* Please connect to app within five minutes after the screen display.

Copyright



Copyright © 2014 FiiSmart Inc. All Rights Reserved.

Fii is the registered trademark of Fiismart.

All Rights Reserved.


All the names of products and companies mentioned in this guide are the trademarks of their respective owners.

Disclaimers


iFever app is a cell phone software developed by FiiSmart Inc., customers use this software on a voluntary basis and should also be aware of the risks. We will not provide any kind of technical support and warranty and will not take any responsibility for using this app unless the customer purchases the iFever device. We will not take any responsibility for information based on this software. Certificate authority based on electronic documents is same as the paper-based ones. They both have the same force of law. By installing 'iFever', you acknowledge that you have read and agreed to the terms of this certification. You are not only granted with the authorities of the above terms but also subject to the restrictions and limitations of it. Any violation of this agreement are directly in breach of the authority and constitute copyright infringement. FiiSmart Inc. has the right to terminate the authority and order to stop the infringements. FiiSmart Inc. reserves the right to investigate and affix the responsibilities.

The correct way of wearing

- Put sensor in the armpit

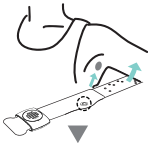

- Clamping

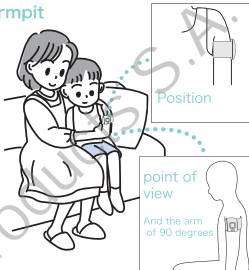

※ For babies please hold the arm, to prevent loosening

- If the position you wear it not correct, might not be able to obtain accurate temperature

Replace the battery



1

Rotate anticlockwise up to alignment bump




2

Open the battery cover



3

Insert the battery

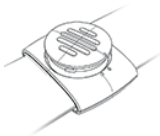

4

Keep and press battery cover and rotate clockwise to alignment bump - to close.

Caring service

We hope that every baby will grow up healthily and we really appreciate every opportunity to serve customers. Please feel free to contact us when you need us.

E-mail: info@fiismart.com

Tel: +86(0755) 86562669

Thermometer is warranted for a period of one year from the date of original retail purchase. Please take one minute to read the following information:

- 1、 Please contact the customer service with your suggestions;
- 2、 Please fill out the service card, which will make communication easier and faster.

The warranty certificate

Date : _____

Add : _____

Name : _____

Supplier: _____

Sign : _____



Scan


FiiSmart Inc.

www.fiismart.com

103, 10 Floor, C1 Building, Nanshan Zhiyuan,
1001 Xueyuan Road, Xili Town, Nanshan District,
Shenzhen, Guangdong Province, China

What is the normal temperature range of body?

Actually all the parts temperature of human body is different, the following form is about the normal temperature range of different body parts.

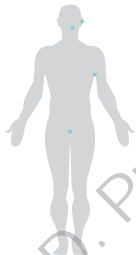

Armpit 34.7°C - 37.3°C

Mouth 35.5°C - 37.5°C

Anus 36.6°C - 38.0°C

Ear 35.8°C - 38.0°C



The correct place to wear iFever thermometer is the upper arm and under armpit, pls be kindly noted that it would affect the accuracy of iFever when wear on another place of body.

QUICK GUIDE

1. Press the button, Wake up device.
2. Download and install the iFever app in App Store or Google Play
3. Register New account by Email.
4. Turn on the Bluetooth function in the smart phone.
5. Select Baby on the iFever app menu.
6. Select connect to iFever
- 7 Your devices and iFever should now be connected



Warm note

1. After connection with iFever and phone, even screen is off, iFever device is still working.
2. For leaving long time or abnormal situation, bluetooth is disconnected, pls connect iFever device and phone via pressing the button and click on app, again.
3. It is available to check the temperature date via registering the same account in another place.